

CORIO CD-601F 1/4

CORIO CD-601F Refrigerated/Heating Circulator


The new generation of laboratory circulators for routine work and standard tasks

Refrigerated / Heating Circulators of the new CORIO[™] series distinguish themselves with a great price-to-performance ratio. They are ideal for all standard tasks and routine work in laboratories and industry.

Your advantages

- · Models for internal and external applications
- Bright, white, easy to read display
- Very quiet
- · Easy pump change-over (internal and external circulation)
- External pump connections (M16x1)
- USB connection
- Space-saving cooling coil design yields more usable space in the bath tank
- Bath lid and drain tap included
- Removable ventilation grid
- Refrigeration unit without side vents
- Class III (FL) according to DIN 12876-1



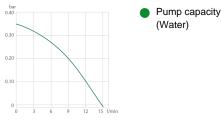

Technical Data

9012705
CORIO
Refrigerated - Heating Circulators
-40 +150
PID
±0.03
0.01 – 0.1 °C
LED
2
°C 20 0 -10 -20 -30 -40 kW 0.6 0.5 0.35 0.2 0.07 0.01
15
0.35
22 x 15 / 20
M16x1
8 / 12
8 10
R452A
150



Global Warming Potential for R452A	2140
Carbon dioxide equivalent stage 1 (t)	0.321
Ambient temperature	540 °C
Dimensions W x L x H (cm)	36 x 46 x 74
Weight (kg)	38.5
Included with each unit	2 portagomma per tubi con d.i. da 8 e 12mm. (Raccordi pompa con filettatura esterna M16x1)
Bath tank	Stainless steel
Bath cover	integrated
Cooling of compressor	Air
Available voltage versions	230 V / 50 Hz 208 V - 230 V / 60 Hz 115 V / 60 Hz 100 V / 50-60 Hz

Characteristics



100% Checked.

100% testing. 100% quality. Each CORIO™ undergoes thorough quality testing before leaving the factory.



ATC.

Absolute Temperature Calibration, 1-point calibration (CD).



Safety.

Quick start.

CORIO CD complies with Class III (FL) according to DIN 12876-1 and switches off automatically in case of high temperature or low liquid level alarm.

Individual JULABO consultation and

comprehensive manuals at your disposal.



Brilliant.

Bright display with strong illumination can be viewed easily from a long distance.



Condensation protection.

Superb design solution. Integrated ventilation directs air over the bath lid and minimizes condensation.



Connection. Easy.

Inclined pump connections (M16x1) facilitate the connection of applications. Each unit includes 2 barbed fittings of 8/12 mm diameter each.




Satisfied customers.

11 subsidiaries and more than 100 partners worldwide guarantee fast and qualified JULABO support.



Services 24/7.

Around the clock availability. You can find suitable accessories, data sheets, manuals, case studies, and more at www.julabo.com.

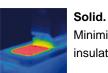

Minimized energy loss through high-quality insulation.



Product data sheet



Everything at the front.

All operating controls and safety functions are accessed easily and comfortably from the front.



Space saving. Free up space.

Place your CORIO right next to an application, another unit, or wall. That saves space. This is made possible by eliminating vents and connections on the sides.

Rubber feet allow for a secured footing of your CORIO[™] to prevent damage to your



Exact.

You can rely on it. PID1 control and 'Active Cooling Control' make the new CORIO™ precise and perfect.



Green technology. Development consistently applied

environmentally friendly materials and technologies.



Handle with ease.

Makes day-to-day work easy. Comfortably move your CORIO[™] around by using the ergonomic handles (front and rear).



Internal. External.

The pump is controlled via a lever located directly below the display. Easily change between internal and external circulation.



JULABO. Quality.

Highest standards of quality for a long product life.



Locked in. The lockable power plug guarantees safe connection. More process safety.



Mobile.

Extra easy handling. Integrated castors for easy repositioning of refrigerated circulators. (CD-600F, CD-601F, CD-900F, CD-1000F, CD-1001F)



Switch on. And off you go. Intelligent operating concept. Ready for operation with just a few quick and easy



Switch on. And off you go. Intelligent operating concept. Ready for

operation with just a few quick and easy steps.



Tidy.

Stable.

steps.

laboratory equipment.

The special drain tap for easy draining of bath fluids without tools.



Timer. Integrated.

CORIO circulators include an integrated timer function. When the set time has elapsed, a signal sounds and the device switches off. Setting range: 0 ... 999 minutes



Touching permitted. Optimum safety. The ergonomic plastic

handle protects your fingers from hot surfaces.

Unique.

Specially designed pyramids ensure that the condensed liquid drips directly back into the bath.





Product data sheet



Modern. Reliable.

High-grade components in every CORIO – platinum sensors, proven motor technology, CAN-Bus communication and much more.



USB.

Remote control made easy. CORIO CD circulators feature a USB connection.

Refer to www.julabo.com for more information regarding the entire JULABO product portfolio. Technical changes without prior notification. Images may deviate from the original.