

CORIO CD-B5 Heating Circulator with Open Bath


The new generation of laboratory circulators for routine work and standard tasks

The Open Heating Bath Circulators of the new CORIO™ series feature high-quality bath tanks made of stainless steel and pump connections.

Your advantages

- · Models for internal and external applications
- Bright, white, easy to read display
- Very quiet
- Easy pump change-over between internal and external circulation
- USB connection
- · High-quality bath tanks made of stainless steel
- Integrated drain screw (except B5)



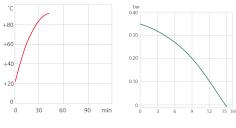

Heat-up time (Water)

Pump capacity (Water)

Technical Data

Order No.	9012405
Model series	CORIO
Category	Heating Circulators with Open Bath
Working temperature range (°C)	+20 +150
Temperature stability (°C)	±0.03
Setting / display resolution	0.01 – 0.1 °C
Temperature Display	LED
Heating capacity (kW)	2
Pump capacity flow rate (I/min)	15
Pump capacity flow pressure (bar)	0.35
Bath opening / bath depth (W x L / D cm)	15 x 15 / 15
Pump connections	M16 x 1
Filling volume liters	3.5 5
Ambient temperature	540 °C
Dimensions W x L x H (cm)	23 x 38 x 41
Weight (kg)	8.2
Included with each unit	2 barbed fittings suitable for tubing 8 and 12 mm inner dia. (pump connections M16x1 male)



Cooling coil	optional
Bath tank	Stainless steel
Bath cover	optional
Available voltage versions	230 V / 50 Hz 208 V - 230 V / 60 Hz 115 V / 60 Hz 100 V / 50-60 Hz

Characteristics



100% Checked.

100% testing. 100% quality. Each CORIO™ undergoes thorough quality testing before leaving the factory.



Modern. Reliable.

High-grade components in every CORIO platinum sensors, proven motor technology, CAN-Bus communication and much more.



ATC.

Absolute Temperature Calibration, 1-point calibration (CD).



Quick start.

Individual JULABO consultation and comprehensive manuals at your disposal.



Brilliant.

Bright display with strong illumination can be viewed easily from a long distance.



Safety.

CORIO CD complies with Class III (FL) according to DIN 12876-1 and switches off automatically in case of high temperature or low liquid level alarm.



Connection. Easy.

Inclined pump connections (M16x1) facilitate the connection of applications. Each unit includes 2 barbed fittings of 8/12 mm diameter each.



Satisfied customers.

11 subsidiaries and more than 100 partners worldwide guarantee fast and qualified JULABO support.



Everything at the front.

All operating controls and safety functions are accessed easily and comfortably from the



Services 24/7.

Around the clock availability. You can find suitable accessories, data sheets, manuals, case studies, and more at www.julabo.com.



Green technology.

Development consistently applied environmentally friendly materials and technologies.



Switch on. And off you go.

Intelligent operating concept. Ready for operation with just a few quick and easy steps.



Internal. External.

The pump is controlled via a lever located directly below the display. Easily change between internal and external circulation.



Timer. Integrated.

CORIO circulators include an integrated timer function. When the set time has elapsed, a signal sounds and the device switches off.



Setting range: 0 ... 999 minutes



JULABO. Quality.

Highest standards of quality for a long product



Unique.

Specially designed pyramids ensure that the condensed liquid drips directly back into the



Locked in.

The lockable power plug guarantees safe connection. More process safety.



USB.

Remote control made easy. CORIO CD circulators feature a USB connection.