# **CORIO CD-BC6 Heating Circulator**


# The new generation of laboratory circulators for routine work and standard tasks

CORIO<sup>™</sup> heating circulators feature professional technology for high standards. These units facilitate internal temperature control in the bath tank or control of externally connected applications.

#### Your advantages

- · Models for internal and external applications
- Bright, white, easy to read display
- Very quiet
- Easy pump change-over between internal and external circulation
- USB connection
- High-quality bath tanks made of stainless steel with bath lid and drain tap
- · Integrated cooling coil for counter-cooling



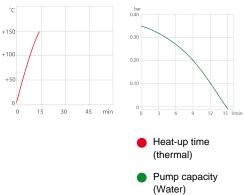

#### **Technical Data**

9012506
CORIO
Heating Circulators
+20 +150
±0.03
0.01 – 0.1 °C
LED
2
15
0.35
13 x 15 / 20
M16x1
8 / 12
4.5 6
540 °C
24 x 44 x 47
10



2 barbed fittings suitable for tubing 8 and 12 mm inner dia. (pump connections M16x1 male)
integrated
Stainless steel
integrated
230 V / 50 Hz 208 V - 230 V / 60 Hz 115 V / 60 Hz 100 V / 50-60 Hz

#### Characteristics



#### 100% Checked.

100% testing. 100% quality. Each CORIO<sup>™</sup> undergoes thorough quality testing before leaving the factory.



# ATC.

Absolute Temperature Calibration, 1-point calibration (CD).



#### Brilliant.

Bright display with strong illumination can be viewed easily from a long distance.



### Connection. Easy.

Inclined pump connections (M16x1) facilitate the connection of applications. Each unit includes 2 barbed fittings of 8/12 mm diameter each.



**Everything at the front.** All operating controls and safety functions are accessed easily and comfortably from the front.



**Green technology.** Development consistently applied environmentally friendly materials and technologies.



# Safety.

CORIO CD complies with Class III (FL) according to DIN 12876-1 and switches off automatically in case of high temperature or low liquid level alarm.



### Satisfied customers.

11 subsidiaries and more than 100 partners worldwide guarantee fast and qualified JULABO support.



#### Services 24/7.

Around the clock availability. You can find suitable accessories, data sheets, manuals, case studies, and more at www.julabo.com.



#### Solid.

Minimized energy loss through high-quality insulation.



#### Stable.

Rubber feet allow for a secured footing of your CORIO<sup>™</sup> to prevent damage to your laboratory equipment.



## Switch on. And off you go.

Intelligent operating concept. Ready for operation with just a few quick and easy steps.



Product data sheet

# CORIO CD-BC6 3/3



Handle with ease.

Makes day-to-day work easy. Comfortably move your CORIO<sup>™</sup> around by using the ergonomic handles (front and rear).



#### Internal. External.

The pump is controlled via a lever located directly below the display. Easily change between internal and external circulation.



# The special drain tap for easy draining of bath fluids without tools.



#### Timer. Integrated.

CORIO circulators include an integrated timer function. When the set time has elapsed, a signal sounds and the device switches off. Setting range: 0 ... 999 minutes



JULABO. Quality. Highest standards of quality for a long product



## Touching permitted.

Optimum safety. The ergonomic plastic handle protects your fingers from hot surfaces.



## Locked in.

life.

The lockable power plug guarantees safe connection. More process safety.



#### Modern. Reliable.

High-grade components in every CORIO – platinum sensors, proven motor technology, CAN-Bus communication and much more.



#### Unique.

Specially designed pyramids ensure that the condensed liquid drips directly back into the bath.



#### USB.

Remote control made easy. CORIO CD circulators feature a USB connection.

Refer to www.julabo.com for more information regarding the entire JULABO product portfolio.

Technical changes without prior notification. Images may deviate from the original. For applications near or below ambient temperature: use a cooling coil or JULABO immersion cooler.