F25-HL Refrigerated/Heating Circulator

Superior models for most demanding applications

JULABO Refrigerated & Heating Circulators for heating and cooling are suitable for external temperature tasks and for temperature applications directly in the bath. A model from three different series' with different capacities can be selected. Powerful circulating pump systems and high heating and cooling capacities guarantee short heat-up and cool-down times. The inevitable dust accumulation that occurs in refrigeration systems is solved by means of an easily removed venting grill, which allows for easy cleaning. Depending on the model, the units are equipped with handles or castors for an easy transportation. All models have a drain tap on the front allowing the fluid to be drained easily. The compact assembly of each unit enables a space to be saved in the laboratory.


Your advantages

- VFD COMFORT DISPLAY
- LCD DIALOG DISPLAY backlit for convenient interactive operation
- · Keypad for setpoints, warning/safety values and menu functions
- ICC (Intelligent Cascade Control), self-optimizing temperature control
- TCF Temperature Control Features to optimize the control behaviour
- ATC3 3-Point-Calibration
- Pt100 External sensor connection for measurement and control
- SMART PUMP, electronically adjustable pump stages
- Druck- und Saugpumpe
- Early warning system for low liquid level (DBGM 203 06 059.8)
- · Adjustable high temperature cut-out, visible via display
- RS232/RS485 interface for online communication
- Optional: Analog interface
- Integrated programmer for 6 x 60 program steps
- · Connections for solenoid valve and HSP booster pump
- Active Cooling Control

Technical Data

Order No.	9312625
Model series	HighTech
Category	Refrigerated - Heating Circulators
Working temperature range (°C)	-28 +200
Temperature control	ICC
Temperature stability (°C)	±0.01
Setting / display resolution	0.01 °C
Integrated programmer	6x60 steps
Temperature Display	VFD, LCD



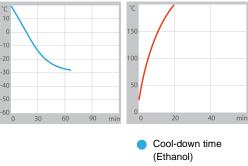





Heating capacity (kW)	2		
Cooling capacity (Medium Ethanol)	°C 20 0 -20 kW 0.26 0.2 0.06		
Pump capacity flow rate (I/min)	22-26		
Pump capacity flow pressure (bar)	0.4-0.7		
Pump capacity flow suction (bar)	0.2-0.4		
Bath opening / bath depth (W x L / D cm)	12 x 14 / 14		
Pump connections	M16x1		
Barbed fittings diameter (inner dia. / mm)	8 / 12		
Filling volume liters	4.5		
Refrigerant stage 1	R134a		
Filling volume refrigerant stage 1 (g)	84		
Global Warming Potential for R134a	1430		
Carbon dioxide equivalent stage 1 (t)	0.12		
External Pt100 sensor connection	integrated		
Digital interface	RS232, RS485 Optional: Profibus		
Ambient temperature	540 °C		
Dimensions W x L x H (cm)	23 x 42 x 64		
Weight (kg)	32		
Classification according to DIN12876-1	Classification III (FL)		
Included with each unit	2 each barbed fittings for tubing 8 and 12 mm inner dia. (pump connections M16x1 female).		
Cooling of compressor	Air		
Power requirement V / Hz / A	230/50/12		
Available voltage versions	230 V / 50 Hz 115 V / 60 Hz 100 V / 50-60 Hz 230 V / 60 Hz		

Characteristics

Display

A perfect view

Ample, easy to read VFD Comfort display for simultaneous display of 3 values, warning functions, high temperature cut-off, pump stages (resolution 0.01 °C)

Additional plain text information Comfortable LCD dialog display for interactive operation with plain text display



Pump stage and liquid level Backlit indicator for selected pump stages and filling volume on Presto® PLUS, Magnum 91 & Forte HT

Operation



Comfortable and detailed Comfortable keypad with additional menu functions for pump stages, calibration, control parameters, programmer, warnings, etc.



Temperature Control

ICC For perfect results

'Intelligent Cascade Control', automatic & self optimizing adjustment of PID control parameters, temperature stability ±0.005 °C internal, <±0.05 °C external

TCF Full control

'Temperature Control Features' for individual optimization, access to all important control parameters, additional settings for band limit, limits, Co-Speedfactor etc.

Highest measuring accuracy 'Absolute Temperature Calibration' for manual compensation of a temperature difference, 3-point calibration

Refrigeration Technology

temperatures

- Easily removable venting grid for quick and easy cleaning
- ACC 100 % Cooling capacity 'Active Cooling Control' for cooling available throughout the entire working temperature range, fast cool-down even at higher


-

Condensation and ice protection A heated cover plate prevents condensation or ice build-up in the bath

Technical Features

Clever pump system

Reliable and consistent pump capacity, electronically adjustable pump stages

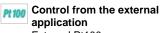

External Pt100 sensor connection for precise measurement and control directly in the external application

Connection compliant to standard RS232/RS485 dual-interface for serial data transmission according to EIA-485 industry standard (2-wire bus technology), upgradable with Profibus DP

//			
1		Þ	
1	υ.		U

Optimal program control For the execution of time and temperature dependant profiles, 6 temperature applies with 60 atoms

temperature profiles with 60 steps max., with real time clock



Connection of additional

equipment Stakei connections for solenoid valve, HSP booster pump and HST booster heater

Warning & Safety Functions



Early warning system for low liquid level Maximum safety for applications, optical and audible alarm, allows user to refill bath fluid before the unit shuts down



Early warning system for high/low temperature limits Maximum safety for applications,

optical and audible alarm, convertible to automated cut-off function



Enhanced protective functions Maximum safety, adjustable high

temperature cut-off or dry-running protection, additional display of setpoints permits easy and precise adjustments



For flammable bath fluids Classification III (FL) according to DIN 12876-1

Refer to www.julabo.com for more information regarding the entire JULABO product portfolio. Technical changes without prior notification. Images may deviate from the original.