F34-HE Refrigerated/Heating Circulator


The powerful upper class unit for demanding applications.

JULABO Refrigerated & Heating Circulators for heating and cooling are suitable for external temperature tasks and for temperature applications directly in the bath. A model from three different series' with different capacities can be selected. Powerful circulating pump systems and high heating and cooling capacities guarantee short heat-up and cool-down times. The inevitable dust accumulation that occurs in refrigeration systems is solved by means of an easily removed venting grid, which allows for easy cleaning. Depending on the model, the units are equipped with handles or castors for an easy transportation. All models have a drain tap on the front allowing the fluid to be drained easily. The compact assembly of each unit enables space to be saved in the laboratory.

Your advantages

- VFD COMFORT DISPLAY
- · Keypad for setpoints, warning/safety values and menu functions
- ICC (Intelligent Cascade Control), self-optimizing temperature control
- ATC3 3-Point-Calibration
- Pt100 External sensor connection for measurement and control
- SMART PUMP, electronically adjustable pump stages
- Adjustable pressure and suction pump
- Early warning system for low liquid level (DBGM 203 06 059.8)
- · Adjustable high temperature cut-out, visible via display
- · RS232 interface for online communication
- · Optional: analogue interfaces
- Integrated programmer for 10 program steps
- Active Cooling Control



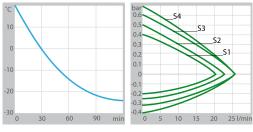

- Cool-down time (thermal)
- Pump capacity (Water)

Technical Data

Order No.	9212634
Model series	HighTech
Category	Refrigerated - Heating Circulators
Working temperature range (°C)	-30 +150
Temperature control	ICC
Temperature stability (°C)	±0.01
Setting / display resolution	0.01 °C
Integrated programmer	1x10 steps
Temperature Display	VFD
Heating capacity (kW)	2
Cooling capacity (Medium Ethanol)	°C 20 0 -20 -30 kW 0.45 0.32 0.14 0.03
Pump capacity flow rate (I/min)	22-26



Pump capacity flow pressure (bar)	0.4-0.7
Pump capacity flow suction (bar)	0.2-0.4
Bath opening / bath depth (W x L / D cm)	24 x 30 / 15
Pump connections	M16x1
Barbed fittings diameter (inner dia. / mm)	8 / 12
Filling volume liters	20
External Pt100 sensor connection	integrated
Digital interface	RS232 Optional: Profibus
Ambient temperature	540 °C
Dimensions W x L x H (cm)	38 x 58 x 64
Weight (kg)	44
Classification according to DIN12876-1	Classification III (FL)
Included with each unit	2 each barbed fittings for tubing 8 and 12 mm inner dia. (pump connections M16x1 female).
Cooling of compressor	Air
Power requirement V / Hz / A	230/50/12
Available voltage versions	230 V / 50 Hz 115 V / 60 Hz 230 V / 60 Hz

Characteristics

Display



A perfect view

Ample, easy to read VFD Comfort display for simultaneous display of 3 values, warning functions, high temperature cut-off, pump stages (resolution 0.01 °C)



Pump stage and liquid level

Backlit indicator for selected pump stages and filling volume on Presto® PLUS, Magnum 91 & Forte HT

Operation



Comfortable and detailed

Comfortable keypad with additional menu functions for pump stages, calibration, control parameters, programmer, warnings, etc.

Temperature Control



For perfect results

'Intelligent Cascade Control', automatic & self optimizing adjustment of PID control parameters, temperature stability ±0.005 °C internal, <±0.05 °C external



TCF Full control

'Temperature Control Features' for individual optimization, access to all important control parameters, additional settings for band limit, limits, Co-Speedfactor etc.



ATC³ Highest measuring accuracy

'Absolute Temperature Calibration' for manual compensation of a temperature difference, 3-point calibration

Refrigeration Technology



Consistent cooling capacity

Easily removable venting grid for quick and easy cleaning



ACC 100 % Cooling capacity

'Active Cooling Control' for cooling available throughout the entire working temperature range, fast cool-down even at higher temperatures

Technical Features



Clever pump system

Reliable and consistent pump capacity, electronically adjustable pump stages



Pt 100 Control from the external application

External Pt100 sensor connection for precise measurement and control directly in the external application



Serial connection

RS232 interface for PC connection, e.g. for data communication and recording of measured values





Easy program control

Integrated programmer for the execution of time and temperature dependant profiles, 1 temperature profile with 10 steps max., with real time clock

Warning & Safety Functions



Early warning system for low liquid level

Maximum safety for applications, optical and audible alarm, allows user to refill bath fluid before the unit shuts down



Early warning system for high/low temperature limits

Maximum safety for applications, optical and audible alarm, convertible to automated cut-off function



Enhanced protective functions

Maximum safety, adjustable high temperature cut-off or dry-running protection, additional display of setpoints permits easy and precise adjustments



For flammable bath fluids

Classification III (FL) according to DIN 12876-1