


ME-31A Visco Bath

Visco Bath for internal temperature applications

JULABO Visco Baths are used for precise temperature applications with capillary viscometers and density meters within the bath tank. The units offer high temperature accuracy and optimal temperature distribution throughout the bath. An integral cooling coil for cooling water is provided for temperature applications below or near ambient temperatures. The units are equipped with a holder for thermometers and therefore offer very good options for temperature calibration procedures and comparisons with accredited thermometers.

Your advantages

- VFD COMFORT DISPLAY
- · Keypad for setpoints, warning/safety values and menu functions
- PID3 cascade temperature control
- ATC3 3-Point-Calibration
- Pt100 External sensor connection for measurement and control
- SMART PUMP, electronically adjustable pump stages
- Early warning system for low liquid level (DBGM 203 06 059.8)
- Adjustable high temperature cut-out, visible via display
- RS232 interface for online communication
- Integrated programmer for 10 program steps

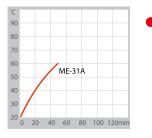

Heat-up tim	ıe
(Water)	

Technical Data

Order No.	9162331
Model series	TopTech
Category	Visco Baths
Working temperature range (°C)	+20 +60
Temperature stability (°C)	±0.01
Setting / display resolution	0.01 °C
Integrated programmer	1x10 steps
Temperature Display	VFD
Heating capacity (kW)	2
Pump capacity flow rate (I/min)	11-16
Pump capacity flow pressure (bar)	0.23-0.45
Bath opening / bath depth (W x L / D cm)	9 x 9 (3 x) / 37
Filling volume liters	31
External Pt100 sensor connection	integrated
Digital interface	RS232 Optional: Profibus
Ambient temperature	540 °C
Dimensions W x L x H (cm)	50 x 20 x 56
Weight (kg)	11
Classification according to DIN12876-1	Classification III (FL)



Cooling coil	integrated
Bath tank	Plexiglass
Number of viscometers	3
Available voltage versions	208 - 230 V / 50 - 60 Hz 100 - 115 V / 50 - 60 Hz

Characteristics

Display



A perfect view

Ample, easy to read VFD Comfort display for simultaneous display of 3 values, warning functions, high temperature cut-off, pump stages (resolution 0.01 °C)



👧 🛅 Pump stage and liquid level

Backlit indicator for selected pump stages and filling volume on Presto® PLUS, Magnum 91 & Forte HT

Operation



Comfortable and detailed

Comfortable keypad with additional menu functions for pump stages, calibration, control parameters, programmer, warnings, etc.

Temperature Control



For higher demands

PID Temperature control with drift compensation and adjustable parameters, improved temperature stability for external applications, temperature stability ±0.01 °C internal, <±0.1 °C external



ATC³ Highest measuring accuracy

'Absolute Temperature Calibration' for manual compensation of a temperature difference, 3-point calibration

Refrigeration Technology



Consistent cooling capacity

Easily removable venting grid for quick and easy cleaning



ACC 100 % Cooling capacity

'Active Cooling Control' for cooling available throughout the entire working temperature range, fast cool-down even at higher temperatures

Technical Features



Clever pump system

Reliable and consistent pump capacity, electronically adjustable pump stages



Pt 100 Control from the external application

External Pt100 sensor connection for precise measurement and control directly in the external application



Serial connection

RS232 interface for PC connection, e.g. for data communication and recording of measured values



Easy program control

Integrated programmer for the execution of time and temperature dependant profiles, 1 temperature profile with 10 steps max., with real time clock

Warning & Safety Functions



Early warning system for low liquid level

Maximum safety for applications, optical and audible alarm, allows user to refill bath fluid before the unit shuts down



Early warning system for high/low temperature limits

Maximum safety for applications, optical and audible alarm, convertible to automated cut-off function



Enhanced protective functions

Maximum safety, adjustable high temperature cut-off or dry-running protection, additional display of setpoints permits easy and precise adjustments



For flammable bath fluids

Classification III (FL) according to DIN 12876-1