

PRESTO A80


temperature control system / process system

Reactor temperature control, tests for all kinds of substances or temperature simulation – the new PRESTO are made for highly precise temperature control and rapid temperature changes. Highly efficient components allow extremely fast compensation of exothermic and endothermic reactions. Lab users benefit from high flow rates, constant pressure, and a controlled build-up of pump pressure. Permanent internal monitoring and self-lubricating pumps contribute to the new PRESTO[®]'s long service life. The integrated 5.7" industrial touch screen displays all important information clearly and concisely enhancing ease of use considerably.

Your advantages

- For highly precise, external temperature applications
- Rapid heating and cooling
- · Wide working temperature ranges without changing fluids
- · Highest performance with small footprint
- · Space-saving design optimizes space utilization in your lab
- NEW 5,7" industrial color TFT touch screen
- NEW USB (Host und Device)
- NEW Ethernet
- NEW SD-Card slot
- RS232 / optional RS485 / optional Profibus DP
- Stand-by input
- · Removable venting grid for simplified removal of dust



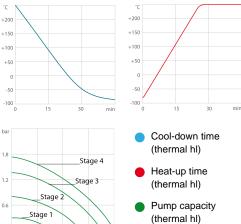

Technical Data

Order No.	9420801
Category	Temperature Control PRESTO
Working temperature range (°C)	-80 +250
Temperature control	ICC
Temperature stability (°C)	±0.01 ±0.05
Setting / display resolution	0.01 °C
Integrated programmer	8x60 steps
Temperature Display	TFT Touchscreen
Heating capacity (kW)	1.8



Cooling capacity (Medium: JULABO Thermal Ethanol)	°C 200 20 0 -20 -40 -60 -80 kW 1.2 1.2 1.1 1.1 0.65 0.1
Pump capacity flow rate (I/min)	16 40
Pump capacity flow pressure (bar)	0.3 1.7
Pump connections	M24x1.5
Refrigerant stage 1	R507
Filling volume refrigerant stage 1 (g)	1140
Global Warming Potential for R507	3985
Carbon dioxide equivalent stage 1 (t)	4.543
Refrigerant stage 2	R23
Filling volume refrigerant stage 2 (g)	500
Global Warming Potential for R23	14800
Carbon dioxide equivalent stage 2 (t)	7.4
External Pt100 sensor connection	integrated
Digital interface	RS232, SD memory card, USB, Ethernet, Modbus Optional: RS485, Profibus
Analog connection input / output	Optional
Ambient temperature	5 40 °C
Dimensions W x L x H (cm)	43 x 65 x 126
Weight (kg)	164
Sound pressure level (distance 1 m) max. (dBA)	68
Process volume min. (active heat exchanger volume) liters	3.9 (1.7)
Internal usable expansion vol. (liters)	5.6
Classification according to DIN12876-1	Classification III (FL)
Cooling of compressor	2-stage Air
Available voltage versions	208V/60Hz (-10/+15%) / 15A / Nema N6-20 Plug 230V/50Hz (+/- 10%) / 13A / UK Plug type BS1363A 230V/50Hz (+/- 10%) / 16A / CEE 7/4 Plug type F

Characteristics

Display


State-of-the-art display technology

TFT Display for comfortable user guidance, colored display of measurement values, graphs and control options, user-defined views

Operation



Optimal ease of use Touch screen for direct operation via display

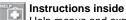


Help menus and explanations in plain text for all control options, help messages and warning messages



Multilingual user guidance

Language selection for display of control options, notifications and warning messages via touchscreen

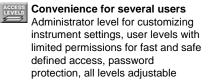



Temperature Control

For perfect results

'Intelligent Cascade Control', automatic & self optimizing adjustment of PID control parameters, temperature stability ±0.01 °C ... <±0.2 °C

TCF Full control

'Temperature Control Features', for individual optimization, access to all important control parameters, additional settings for band limit, limits, co-speedfactor etc.

Control from the external application

External Pt100 sensor connection for precise measurement and control directly in the external application

Highest measuring accuracy ATC³

'Absolute Temperature Calibration' for manual compensation of a temperature difference, 3-point calibration

Refrigeration Technology

6*

Consistent cooling capacity Easily removable venting grid for

quick and easy cleaning

ACC 100 % Cooling capacity

'Active Cooling Control' for cooling available throughout the entire working temperature range, fast cool-down even at higher temperatures

Energy saving cooling

Proportional cooling control for automatic adjustment of cooling power or temporary switch-off of compressor as needed to save up to 90 % energy in comparison to unregulated cooling machines

Technical Features

Intelligent pump system


Reliable and consistent pump capacity, electronically adjustable pump stages or pressure value, automatic adjustment of pump capacity to viscosity

Communication via networks For the remote control of instruments via Ethernet networks, full access to all functions of the unit via a networkcapable PC



Intelligent communication USB connection for data exchange (e.g. service data) or for wireless remote control via WirelessTEMP®

Data exchange via SD-Card For data exchange (e.g. service data) via SD memory card

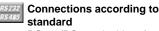

RS232/RS485 dual-interface for serial data transmission according to EIA-485 industry standard (2-wire bus technology), upgradable with Profibus DP

Comfortable program control Integrated programmer for the

execution of time and temperature dependant profiles, 8 temperature profiles with 60 steps max., with real time clock



Quiet as a whisper

Efficient components produce only a minimal sound decibel level





All connections as well supply and exhaust air are located at the front or rear, no venting grids on the sides, units can be placed close to each other or the application



Continuous operation up to +40 °C


Robust temperature control instrument, continuous operation even at ambient temperatures of up to +40 °C



)Em

Easy transport by one person Ergonomic design facilitates moving and positioning by one person

Filling level at a glance

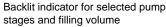

Warning & Safety Functions



Early warning system for high/low temperature limits Maximum safety for applications, optical and audible signal when limits are exceeded.


00 Duplicate safety

Adjustable high temperature cut-off for internal tank and for integrated expansion vessel



For flammable bath fluid

Classification III (FL) according to DIN 12876-1

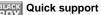

If an error occurs, the integrated Black-Box function permits fast diagnosis by the JULABO service team